



**GLACIER**

**-86°C ULTRALOW TEMPERATURE FREEZERS**



## Environmentally Safe

- HFC (CFC-Free) Refrigerants
- (CFC-Free) Polyurethane Foamed In-Place Insulation
- Low Noise Operation 47 dB(A)
- Efficient Design Compressors



For more information contact NuAire, Inc. at 1.800.328.3352 or [www.nuaire.com](http://www.nuaire.com)


## Controls


- Alarm System
- High / Low Temperature Alarm ( $\pm 5^{\circ}\text{C}$  to  $\pm 20^{\circ}\text{C}$  adjustable)
- Remote Alarm Dry Contact
- Continuous Self Diagnostics (sensors abnormality)
- Battery Check (battery replacement recommended)
- Self Diagnostic Features
- Microprocessor Controls
- Automated  $\text{CO}_2$  Back-up System
- Filter Check (filter replacement recommended)
- Compressor Time Delay Adjustment 2-15 min.
- Temperature Zero Adjustment
- Remote Alarm Delay Adjustment



## Ergonomics and Safety

- Eye Level Control Center
- Single Cam-Action Outer Door Catch
- "Claw" Inner Door Latches **A**
- Halo Bypass Coil Frame Heater (minimizes gasket frost buildup)
- Adjustable Chamber Shelves **B**
- Deadbolt lock with key **C**
- Slide-in-Place Outer Door **D** (lift up feature)
- Audible and Visual Alarms
- Backup Redundant Systems
- Internal Chamber Pressure Release Gasket
- Exterior Door Latch **E** (durable, reliable mechanism)

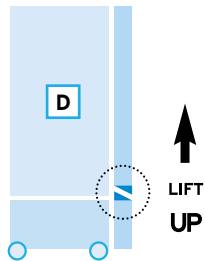

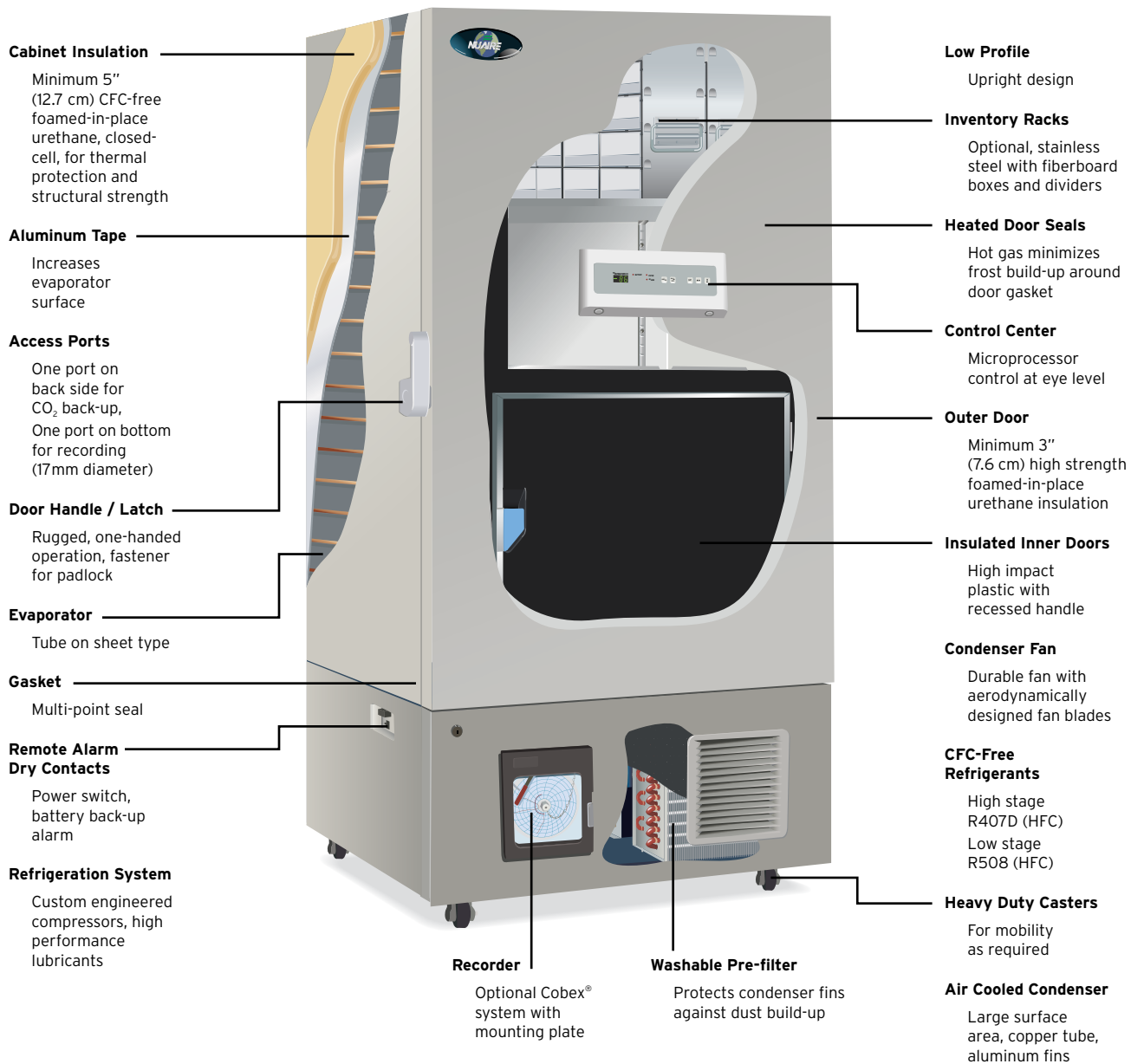

**NU-9000 Series** -86°C Ultralow Temperature Freezers

| Specifications | Model     | Ft. <sup>3</sup> /<br>Liters | Electrical | Dimensions<br>(Exterior W x D x H) | Dimensions<br>(Interior W x D x H) | Shipping<br>Dims.     |
|----------------|-----------|------------------------------|------------|------------------------------------|------------------------------------|-----------------------|
|                | NU-9483GC | 17.1 / 483                   | 220V 60 Hz | 35.0 x 34.5 x 78.3 in.             | 24.8 x 23.6 x 50.4 in.             | 85.8 ft <sup>3</sup>  |
|                | NU-9483E  | 17.1 / 483                   | 230V 50 Hz | 890 x 876 x 1990 mm                | 630 x 600 x 1280 mm                | 2.43 m <sup>3</sup>   |
|                | NU-9668GC | 23.6 / 668                   | 220V 60 Hz | 44.5 x 34.5 x 78.3 in.             | 34.3 x 23.6 x 50.4 in.             | 105.9 ft <sup>3</sup> |
|                | NU-9668E  | 23.6 / 668                   | 230V 50 Hz | 1130 x 876 x 1990 mm               | 870 x 600 x 1280 mm                | 3.00 m <sup>3</sup>   |

\*Specify model with appropriate letter suffix for electrical specifications. "NU-9334E" for 230 VAC / 60 Hz

**Features**

**Standard Features**

- GC model includes voltage booster
- CFC-free foam insulation  
127 mm [5 Inches]
- Access port
- Heavy-gauge steel cabinet
- Heavy-duty swivel casters
- Claw-like inner door latches
- Key lock
- Two [2] air insulated inner doors
- Four [4] compartments
- Four [4] shelves – Three [3] adjustable
- Multi-point gasket seals
- Eye-level controls
- Four [4] year parts only warranty (International)
- Two [2] year parts and labor warranty (U.S. and Canada)

**Optional Features**

- Cobex® temperature recorder
- CO<sub>2</sub> back-up system [NU-UB2]
- Painted steel interior
- Surge suppressor
- Cryo gloves
- Floor anchor system
- Inventory system

**Electrical Requirements**

- 220V 60Hz; 230V 50Hz
- All single phase

**Temperature Control System**

- Microprocessor control system
- Digital temperature display
- Set point security system

**Refrigeration System**

- Air-cooled cascade refrigeration system
- Two [2] air-cooled compressors
- Temperature range -50°C to -86°C
- 100% HFC and CFC-free refrigerants
- Efficient down feed evaporator
- High capacity air-cooled condenser

**Alarm Systems**

- Microprocessor alarm
- Battery back-up
- Remote alarm contacts



CE All models are designed to meet CE and SGS USTC criteria for safety and performance.



## Options and Accessories



### CO<sub>2</sub> Back-Up System

Will maintain -70°C at a 45 Kg. / Hr. rate.



### Cryo Safety Gloves

Protection for Hands and Arms



### Temperature Recorder

Easily Installed into the Freezer Control Panel




### 2" or 3" Fiberboard Boxes

## Freezer Inventory Racks

Freezer racks come pre-packaged with 20 x 2" (51 mm) cardboard boxes with 81 cell dividers per box.


### NU-9483 16 Racks

320 (4 deep, 5 high) 2" Boxes | Part# SV-HL-20-2-81

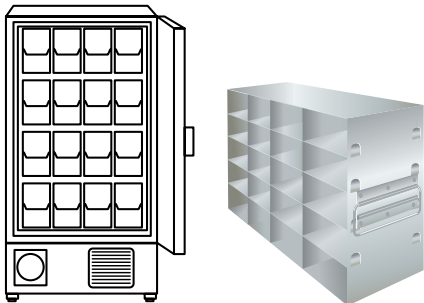

### NU-9668 24 Racks

480 (4 deep, 5 high) 2" Boxes | Part# SV-HL-20-2-81

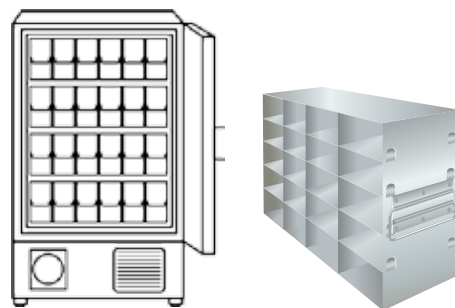

### To Locate Your Local NuAire Sales/Service Representative

Please visit NuAire's Web site at: [www.nuaire.com](http://www.nuaire.com) and select the 'Sales Team' menu item. Information can be found there for both United States and international distributors.



### Other NuAire Products



CO<sub>2</sub> Incubators



Biological Safety Cabinets

Benchtop and Micro Centrifuges



Laminar Airflow Products



For more information please visit  
[www.nuaire.com](http://www.nuaire.com) or call **1.800.328.3352**

