



### Vertical sliding door sterilizers designed for ease of use, simplified maintenance, and low cost of ownership.

The Consolidated EZ-Glide™ Sliding Door Series offers high performance sterilization capacity in a minimal footprint. The sleek vertical door design includes ergonomic, safety, and engineering features that simplify operation, improve serviceability, and document sterilization function for record keeping.

#### EZ-Glide™ Sliding Door




- The EZ-Glide™ sliding door steam sterilizer simplifies chamber access with a sliding counterbalanced door that seals safely against an air actuated non-stick gasket. The non-stick gasket design provides a longer service life than a conventional steam-actuated gasket.
- The gasket is recessed into a precision -machined slot and will maintain a secure seal even in the event of a total utility loss. The gasket is easy to replace and does not require the use of tools for replacement or lubrication for operation.
- The door sealing mechanism includes redundancy for added security; a safety prevents the door from opening if the chamber is under pressure.
- Ergonomic benefits include a wide-grip handle and smooth, vertical open/close action on a counterbalanced assembly that is easy-to-operate with one-hand access.
- Available in both manual and powered door operation
- One touch open feature
- Safety features:
  1. Low closing force, can be stopped with minimal pressure
  2. Powered door is closed using a push and hold button

#### ProGEN™ Steam Generator

When house steam is not available or inadequate, the Consolidated ProGEN™ electric steam generator serves as the steam source. Developed to fit directly below the sterilizer chamber within the chassis, the ProGEN™ is designed to deliver a high rate of steam to the chamber without lag time between cycles. This increased steam capacity accelerates recovery and achieves cycle setpoints more quickly to yield the most efficient turn-around time possible.

The ProGEN™ includes front-access to open-source heating elements for quick replacement over the life of the sterilizer, as well as:

- Larger than normal steam capacity, resulting in faster up-time.
- Large service opening for easy generator maintenance and inspection.
- Whisper-quiet automatic feed water pump, pressure regulator and electronic water level control.
- Standard voltages 208, 240, 380v, 480 in single phase or three phase.
- Capacities from 20 kW to 45 kW.
- Optional automatic blowdown.
- Available for retrofit to existing sterilizer systems.


Model	Chamber Dimensions (w x h x f-b)	Volume (cu. ft. / liters)	Ship Weight (lbs/kg)	Overall Dimensions (w x h x f-b)
Sliding Door Model				
 3AV	 26AV	 26BV		
3AV	20" x 20" x 38" 50.8 x 50.8 x 96.5 cm	8.8 cu. ft. 250 liters	1,650 lbs 750 kg	39.375" x 71" x 47" 100 x 180.3 x 119.4 cm
26AV	26" x 26" x 39" 66 x 66 x 99 cm	15.25 cu. ft. 430 liters	2600 lbs 1182 kg	46.375" x 77.25" x 54" 118 x 196 x 138 cm
26BV	26" x 26" x 49" 66 x 66 x 124.5 cm	19.2 cu. ft. 544 liters	2850 lbs 1296 kg	46.375" x 77.25" x 64" 118 x 196 x 163 cm

Specify Advantage Series™ ADV-PRO™, ADV-PLUS™ or ADV-PB™ when ordering.  
Available with optional ProGEN™ integral electric steam generator.



Model 3AV

## Designed for Maximum Up-Time

Feature	Benefit
Unique non-stick door gasket*	Specially formulated to eliminate sticking and to avoid door jamming
Air-actuated gasket	Minimizes steam exposure to greatly extend gasket working life
Whisper-quiet, integrated air compressor	Minimizes noise, eliminates need for facility air supply, supplies constant air pressure to gasket
Side-mounted controls	Avoids steam exposure to electronics and extends touchscreen life
Front access swing-out assembly	Simplifies serviceability, assures access to all components
Counterbalanced sliding door	Simplifies working footprint, lifts without the complexity of pneumatics or service-prone hydraulics
Low-mounted counter pulleys	Prevents steam exposure to pulley system, minimizes corrosion on pulley cables
One style solenoid valve used throughout sterilizer*	Simplifies maintenance, non-proprietary valve reduces cost
Integrated post-vac drying cycle	Vacuum drying phase eliminates steam and removes moisture from the load (i.e. reduces wet-packs)
Front-mounted chamber drain	Simplifies access to the chamber drain strainer for easy cleaning
Pressure gauge calibration	Permits periodic calibration to maintain accuracy, extend component life
Customizable design*	All sliding door sterilizers are custom made to exact configuration, functionality and options
Slim profile	Fits through standard doorways
Automatic power door operation	Powered door operation allows 1 finger open/close operation. The door can be manually operated in the event of power loss.
 WaterEco® Conservation System**	Significantly reduces sterilizer water usage




\*CSS Exclusive

\*\*Optional

## Advantage Series™ PLC Controls




The EZ-Glide™ can be ordered with any of the Consolidated Advantage-Series Sterilizer Controllers including ADV-PRO™, ADV-PLUS™ and ADV-PB™.

Controller Models	Display	Key Features
ADV-PRO™		6" or 10" Color Touchscreen 22+ Cycles Pre-Vac Capability Password Protection BSL-3 Compatible
ADV-PLUS™		6" Monochrome Touchscreen 12 Cycles Pre-Vac Capability
ADV-PB™		3" Monochrome Touchscreen 2 Cycles

## SteriNET® Remote Access Features

Advantage™ controllers can integrate with a suite of SteriNET® state-of-the-art monitoring, troubleshooting, supervisory control and data acquisition solutions required to meet specific laboratory and/or facility requirements.

- SteriNET® **DataPort**, data downloading to a local computer or facility monitoring system.
- SteriNET® **Connex** for live factory diagnostics.
- SteriNET® **APPs**, third party monitoring applications for mobile devices.
- SteriNET® **SCADA**, connects to a Supervisory Control and Data Acquisition System.

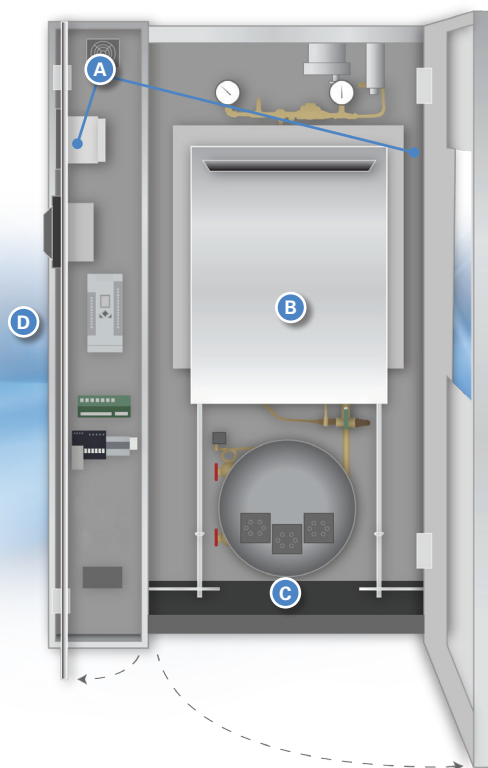

- A.** Hinged front panels for easy access
- B.** EZ-Glide™ vertical sliding door
- C.** ProGEN™ electric steam generator (optional)
- D.** Advantage Series™ control column

**Model 3AV (Front View)** Service access panels on left- and right-hand hinges place all interior components and plumbing within easy reach and front access. Controls, printers, communications, electrical, mechanical and plumbing circuits are positioned for proactive maintenance.



P.O. BOX 297 | 76 ASHFORD ST. | BOSTON, MA 02134-0003

P: 617.782.6072 | F: 617.787.5865 | INFO@CONSTERIL.COM | CONSTERIL.COM

Featuring Green Technology  for Energy Savings and Minimal Environmental Impact



©2014 Consolidated Machine Corporation V1 9/14