



# Mito FB Combi

the comprehensive approach  
for the sanitization of feeding bottles

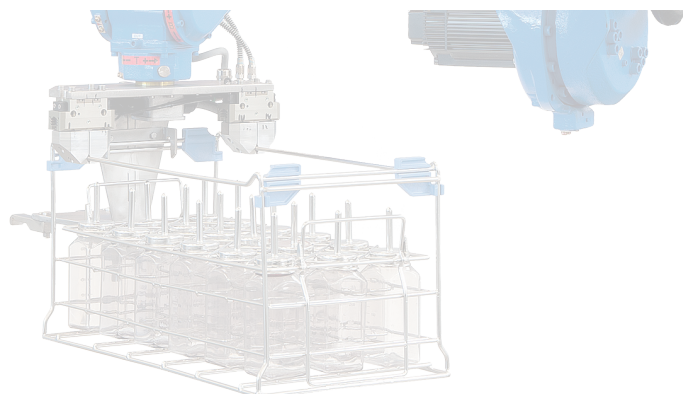

*The most compact, multifunctional and automatized system for the disinfection of feeding bottles and caps at the same time.*

**Decapping, emptying, washing, sanitization, filling and recapping: all in one solution within a small footprint.**

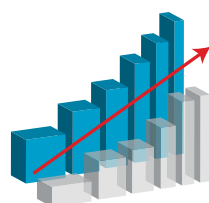



**High productivity  
within a compact  
solution**



- Less than 4 m<sup>2</sup> footprint
- Over 1080 bottles and caps reprocessed per hour with a considerable saving of time and personnel
- Suitable for in line and conventional barrier pass-through applications

## Washer system dimensions

	W	D	H	H+
mm	3400	1150	875	2450
inches	133,8	45,3	35,5	96,5

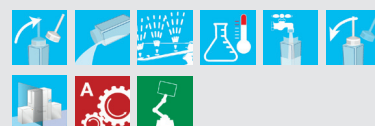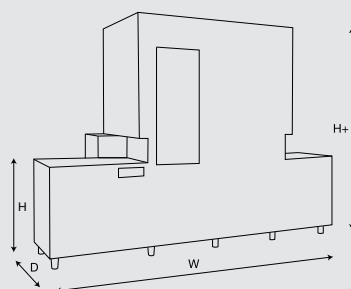

The all in one solution




## MITO FB COMBI

### Automated feeding bottles washing system

**MITO FB COMBI** offers a compact and full solution to support the increasing demand of automated systems for the disinfection of feeding bottles and caps at the same time.

#### Key features

- **High productivity.** Over 1080 bottles processed per hour.
- **Fully automated.** Compatible with robotic loading and unloading functions.
- **Compact solution.** Less than 4 m<sup>2</sup> footprint.
- **Versatile.** MITO FB COMBI is compatible with any model of caps and feeding bottles available on the market.
- Fully automatic process, 1x18 bottle crate capacity.
- PLC control system with 5.7" touch screen LCD display. 3 password level protection programming.
- Chambers entirely made of stainless steel AISI 316L, external body made of AISI 304.
- Full transparent external doors for visual inspection of the process.
- Independent washing and rinsing systems to avoid cross contamination. Individual nozzles are dedicated to outer and inner surface of the bottles.
- Disinfection by rinsing up to 90°C / 194°F.
- Washing and rinsing chambers provided with stainless steel interlocked sliding doors.
- Double HST glass door for washing chamber visual inspection.

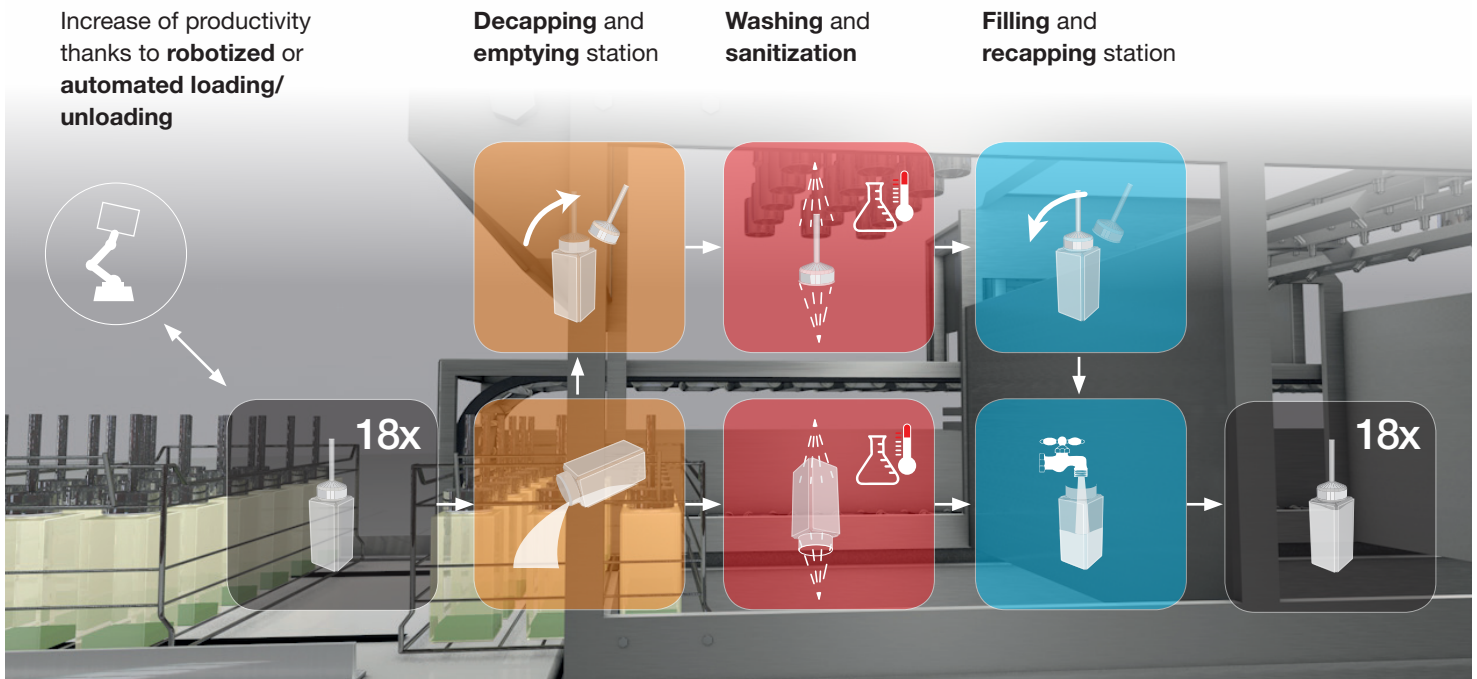

#### Headquarters

##### STEELCO S.p.A.

Via Balegante, 27 - 31039 Riese Pio X (TV) - Italy  
Ph. +39 0423 7561 - Fax +39 0423 755528  
e-mail: steelco@steelcospa.com

[www.steelcospa.com](http://www.steelcospa.com)

#### Branches

##### STEELCO HUNGARY

Temesvár u., 19-21  
1116 Budapest - Hungary  
Ph. +36 1 371 0257  
Fax +36 1 371 0258  
e-mail: steelco@t-online.hu

##### STEELCO USA

999 Stinson Way - Unit 307  
West Palm Beach, FL 33411 USA  
Ph. +1 561 791 8313  
Fax +1 561 791 8213  
e-mail: contact@steelco-usa.com

##### STEELCO ASIA

6th floor, Suite 16, IOI Business Park  
Persiaran Puchong Jaya Selatan - Bandar Puchong Jaya  
47100 Puchong, Selangor - Malaysia  
Ph +60 3 8064 4339 (Ext. 112) - Fax: +60 3 8064 4340  
e-mail: adminasia@steelcospa.com



"For the Environmentally conscious"

